



RAPID INDUCTION ACCURATE MEASUREMENT

- Highly sensitive infrared probe ·
- Noncontact measurement ·
- 1 second can be measured ·

Medical Infrared Thermometer

▶ **AET-R1B1**

LED Digital Display

The black display box focuses the reading, the screen no longer changes color, only bright font. reading clear, and easy to get the temperature value.



1 second measurement

High sensitive probe, the temperature can be measured in 1 second, after a measurement without other operations, you can enter the next measurement it is suitable for multi-person measurement or multi-time measurement.



Authoritative Certification of CE/FDA













Ergonomic Design Comfortable to Hold

Following ergonomic design rules,
the body is solid and stable, handiness,
comfortable grip, firm grasp





Features

| | | | | |
|--|--|---|--|--|
|  a key measure |  memory storage |  replaceable battery |  a second display |  c/f switch |
|  forehead measurement |  hot tip |  LED display |  automatic shutdown |  intelligent accurate |

Technical parameters

Measuring range: 32.0°C~42.2°C

Accuracy: $\pm 0.2^{\circ}\text{C}$ (35.0°C~42.0°C) $\pm 0.4^{\circ}\text{F}$ (95.0°F~107.6°F)
 $\pm 0.3^{\circ}\text{C}$ ($\pm 0.5^{\circ}\text{F}$) other range

Resolution: 0.1°C

Memory: 32

Power supply : 2*AAA(DC 3V)

Dimension: 144mm*101mm*26mm


Weight: 88g(exclude battery)

Auto power-off: 60s

PACKING SIZE

| | |
|--------------|---------------------|
| ■ Qty./CTN | 80pcs |
| ■ gift box | 40*90*152mm |
| ■ carton box | 430*385*330mm |
| ■ CBM | 0.055m ³ |
| ■ G.W | 9KG |

Packaging Accessories

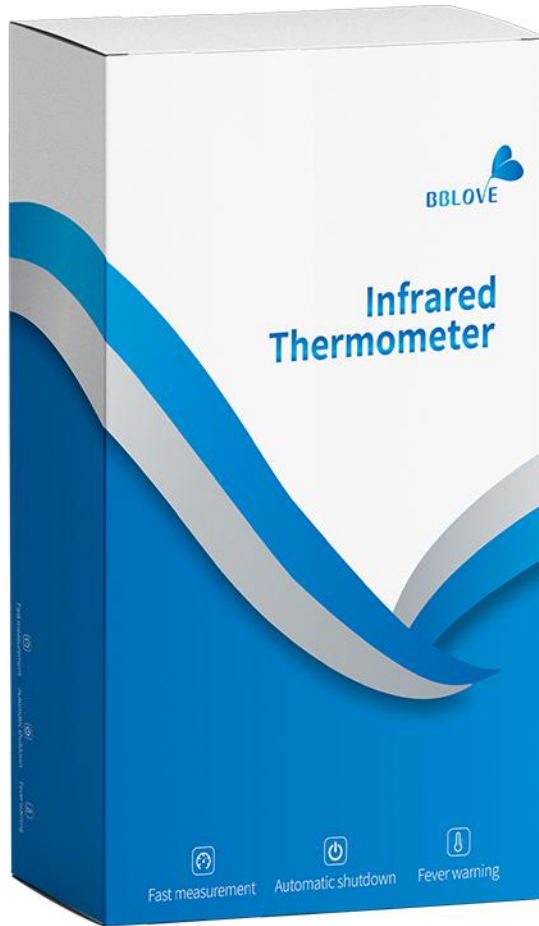

合格证 CERTIFICATE

产品名称:红外额温计
Product name:Infrared Thermometer
型号/规格(Model / Specification):AET-R1B1
执行标准(Executive standard):
YY0505-2012, GB9706.1-2007, GB/T 21417.1-2008
生产许可证编号(Production license number):
粤食药监械生产许20162752号
注册证编号/产品技术要求编号(Registration certificate No. /
product technical requirement No.):粤械注准20162200149

生产日期(Production Date):
生产批号(Production batch No.):

检验员(Inspector):
检验结果:合格
Inspection result:Qualified
使用期限:5年
Service life:5 years
厂家:深圳市爱立康医疗股份有限公司
Manufacturer:Alicon Medical Shenzhen, Inc.
地址:深圳市坪山区坪山街道锦田路14号深福保现代光学产业园B栋第4层
Address:4th Floor, Building B, Shenfubao Modern Optics
factory, No.14, Jinxiu Middle Road, Kengzi Street, Pingshan
District, Shenzhen

Customer retention (please show this card when repairing)

Customer Name: _____
Address: _____
Tel/HP: _____
Date of produce: _____
Model: _____
Batch NO.: _____
Date of purchase: _____

Please keep this warranty card properly. It will not be reissued if lost.



Intended use

Thank you for purchasing our Infrared Thermometer, type AET-R1B1. This Infrared Thermometer is a non-sterile, reusable, handheld device. It can be used by consumers in homecare environment and doctors in clinic as reference. It is intended for measuring human body temperature of three months and above of people by detecting infrared heat from the center of the forehead.
-For safe and proper use of this product, be sure to read and fully understand Safety Precautions contained in this user manual.
-Keep this manual with you for a convenient reference.
-Keep this manual in a proper way and avoid any loss.
-If you need other information, please contact the manufacturer (See the PRODUCT INFORMATION).



Quick Start Guide



The correct measurement method to follow:
Measure the forehead.
Distance 15mm to 50mm

Correct measurement method

1. Turn the thermometer on by pressing the "SCAN/ON" button. To choose mode press "SCAN/ON" on the shutdown mode. A beep sound will be heard following which the screen will indicate with "°C" indicating the unit is ready for use.
2. Place the thermometer probe point near the forehead and press "SCAN/ON" and release after ~ 1 second and the thermometer will beep and display the accurate temperature. Whenever consecutive readings are required, please wait for at least 5 second. The above procedure can be repeated as many times as required.
3. To turn off the device press and hold the "SCAN/ON" switch for ~ 5 seconds and the thermometer will turn off automatically after ~ 60 seconds.



ALIEN

Alicn medical shenzhen Inc.

www.alicn-med.com